



**Curved Small Head to Cut
Rectum in Deep Pelvic
Rapid, Reliable, Viewable**



Short in width, long in curved cutting line

Single Patient Use Transverse Cutter and Reload

Major Features

Great Operation Experience

- Curved small head, easy to enter deep pelvic under direct vision
- Intermediate position, easy to adjust cutting and stapling location
- Separate closing and firing handles, clear operation processes

Great Efficiency

- Short reload width, long cutting line
- Simultaneous cutting and stapling, simplified procedures
- Quick and convenient operation

Great Reliability

- Stapling line 4.5mm longer than cutting line
- Tactile and audible feedback in operation process
- Positioning pin, ensured stapling of all target tissue
- Parallel closure, identical formation of all staples


Golden Stapler Surgical Co., Ltd.



**Advancing Surgical Technology
to Benefit More People**



Curved Small Head, 62mm Width

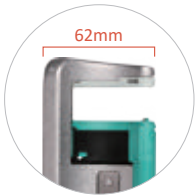

Easier for vertical placement into deep pelvic
Procedure under surgeon's direct vision
4 rows of staples released during firing,
tissue in the middle to be cut
Stapling while cutting, easy and convenient
Reduce surgery processes, more effective

Positioning Pin




Close and lock target tissue, prevent
tissue slippage
Assist parallel closure
Uniform B-shape formation of all staples
Great stapling and hemostasis effect

Curved Reload



Short in width, long in curved cutting line
Compress tissue up to 40mm
Stapling line 4.5mm longer than cutting line
Great stapling and hemostasis in cutting
line's end

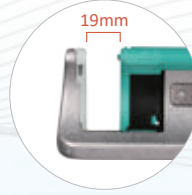

Jaw Opening 19mm Wide


Big opening, easy to put tissue inside

Separate Closing and Firing Handles



Independent closing and firing processes
Prevent miss-operation

Intermediate Position

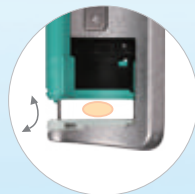

Easy to move the anvil to deep pelvic
Easier to adjust cutting and stapling location



Cutting Washer

Tactile and audible feedback during firing
Clear suggestion of operation process

Code	Color	Cutting Line	Open Staple Height	Staple Rows	Closed Staple Height	Staple Qty	Tissue Thickness
CCSA40G / JYH(B)-40H	Green ●	40mm	4.5mm	4	2.0mm	46	Thick

* Staple from titanium material

Golden Stapler Surgical Co., Ltd.

Add: Building 7A, Jiangsu Wujin Sci-Tech Pioneer Park, 256 Mid Mingxin Road,
Hutang Town, Wujin District, Changzhou, Jiangsu Province, China 213164
Tel: +86 519 85070251 Fax: +86 519 85070252
Web: <http://www.goldenstapler.com>