



Smaller Anvil Head Vertically Makes Access to Lower Pelvic or Lung Hilum Easier and Safer


Single Patient Use Transverse Cutter and Reload

- o Transecting Tissues after Bilaterally Stapling To Make Anastomosis Reliable
- o Smaller Anvil Head Vertically Makes Access to Lower Pelvic Easier
- o Suitable for Colorectal, Gastrointestinal and Thoracic Surgeries
- o Breakthrough Cutting Washer Provides Audible and Tactile Feedback on Every Firing
- o Parallel Jaw Closure Provides Consistent Staple Formation
- o One Handed Operation
- o Four Different Cartridges Reloadable To The Device

Golden Stapler Surgical Co., Ltd.



*Advancing Surgical Technology
to Benefit More People*

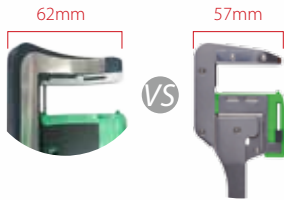

Smaller Anvil Head

Vertically Lower Pelvic Access
Make LAR or Ultralow LAR Easier and Reliable



Breakthrough Cutting Washer

Help Cut Rectum Off
Help Cut Bronchus or Duodenum Off
Deliver Audible And Tactile Feedback
Indicating The Firing Sequence Is Complete



Large Jaw Aperture

Help Obtain Placement Below Tumor Site
Facilitates Maneuverability and Mobility
along the Rectum



Retaining Pin

Captures Intended Tissue before closing device
Intermediate Locking Position Allow Repositioning of the Tissue
Parallel Jaw Closure Provides Consistent Staple Formation



Wave Spring Mechanism

Staples Driven Out First, Knife Out Second
Stapling First, Cutting Second



Closure and Firing Trigger

Close First, And Then Fire
One Handed Operation

PRODUCT CODE	CARTRIDGE CODE	CARTRIDGE COLOR	NUMBER OF STAPLED ROWS	TISSUE THICKNESS	OPEN STAPLE HEIGHT	CLOSED STAPLE HEIGHT	EFFECTIVE STAPLES LINE
CCZB40G / JYH(A)-40H	CCZBR40G	GREEN ●	4	THICK	4.5	2.0mm	40mm

* Staple from titanium material